


Custom Ultrasonics, Inc.

Sediment and Bacteria Water Filtration Assembly
For ONE HOT water line -- (Nominal temp= 110 degF +/- 10)
Two set-ups are required for Hot and Cold water attachment.

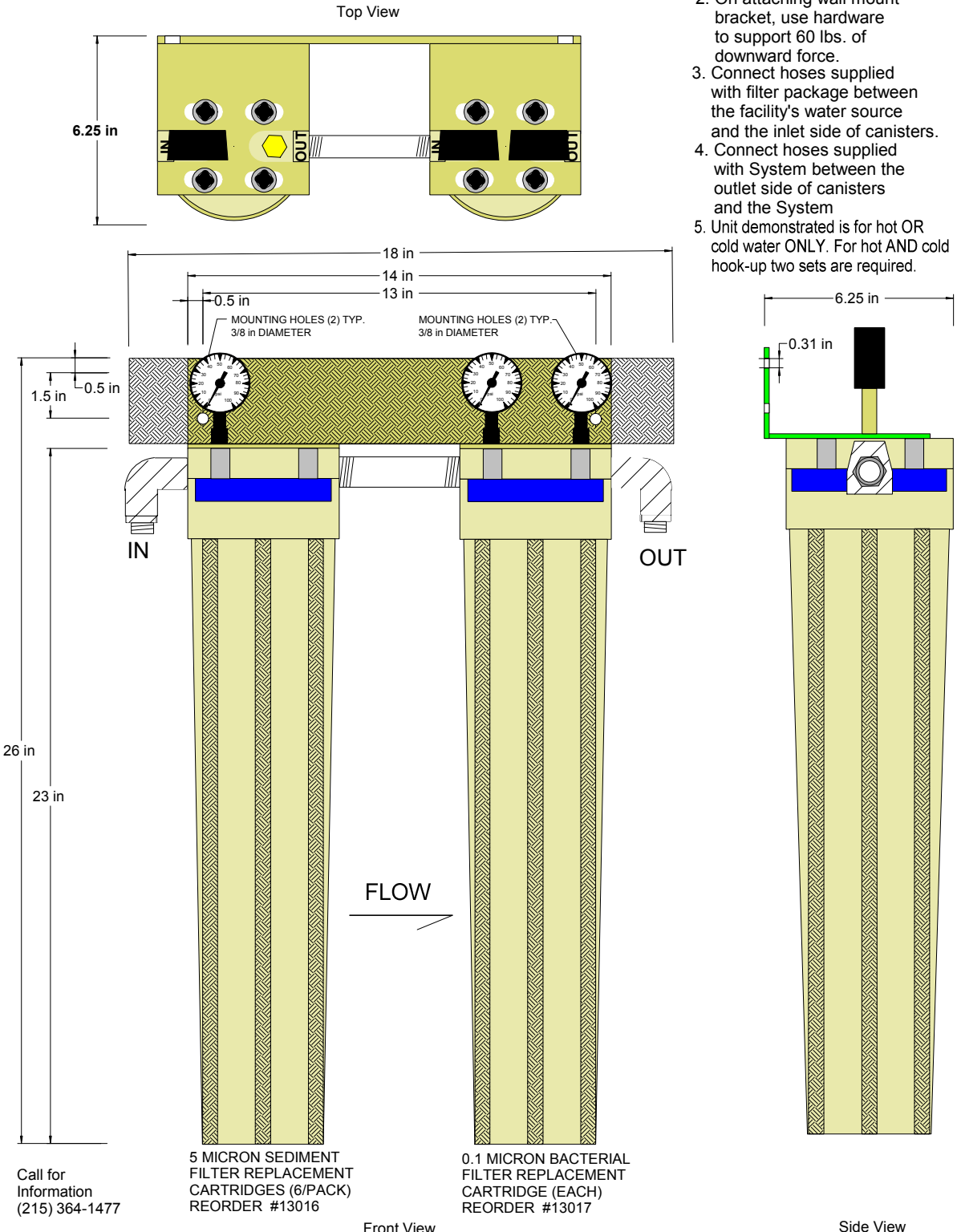

Installation Instructions:

1. Locate the filter housings behind the System approx. 44" AFF [nominal]
2. On attaching wall mount bracket, use hardware to support 60 lbs. of downward force.
3. Connect hoses supplied with filter package between the facility's water source and the inlet side of canisters.
4. Connect hoses supplied with System between the outlet side of canisters and the System
5. Unit demonstrated is for hot OR cold water ONLY. For hot AND cold hook-up two sets are required.

Call for Information
 (215) 364-1477